




User Manual

MYNT DREAMYNT EYE MASSAGER

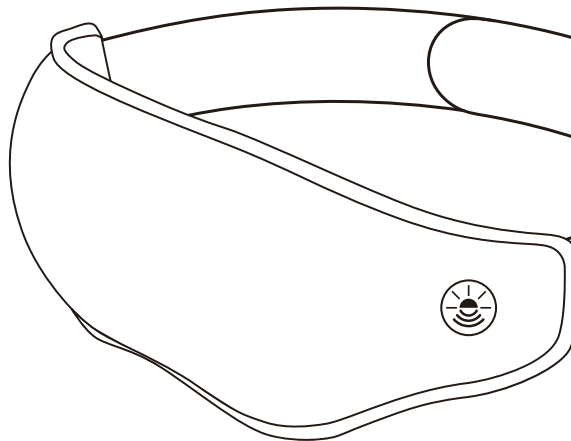

MN-EM02.

Please read and keep this manual carefully before using the product.

ABOUT

The ultra-soft and light-weighted eye mask fits perfectly over your eyes. It provides warm compression as well as high vibrating massage, helping to accelerate the blood circulation. Perfect for eye muscle fatigue.




Using with gel ice pack, it will help to reduce appearance of puffiness and relax the eyes.

Foldable & portable, it's wireless for use anywhere and anytime when you desire an eye massage.


SPECIFICATION

Product Name	MYNT DREAMYNT EYE MASSAGER
Model	MN-EM02
Operating Voltage	3.7V=
Working Power	6W
Charging Voltage	5V==
Charging Current	1A
Heat Temperature	45 ± 5°C
Battery Capacity	600mAh
Product Dimension	L240 x W75 x H95 mm
Net Weight	77.5g

PARTS IDENTIFICATION:

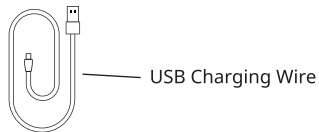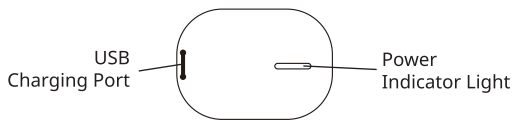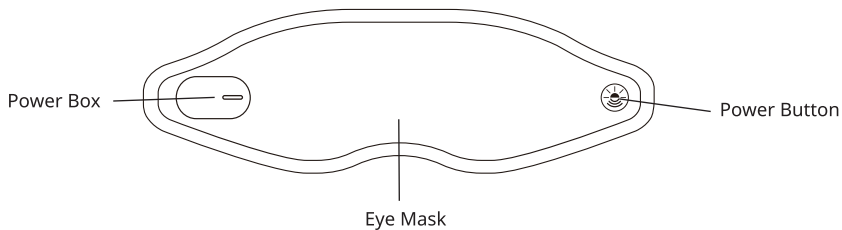

HOW TO USE:

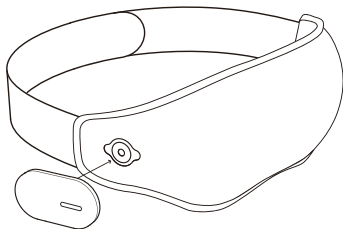


1


Put the power box magnetically attached to the eye mask. After correct installation, the power indicator light will turn on.

2

Fit the eye mask over your eyes, then press and hold the power button on the left side of the eye mask for 1 second. The eye mask will be turn on and there will be a prolonged vibration to indicate.

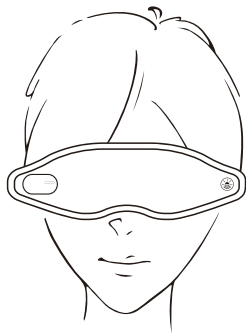

HOW TO USE:

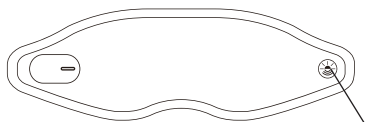


Power Button

3

Power on the default massage mode (heat + vibration), short press the power button to switch between massage mode (heat + vibration) and heat mode (heat only), there will be 2 short vibrations to indicate.

4

If you need cold compress, please refrigerate the gel ice pack for 10 minutes, and then place it inside the eye mask for use without turning on the eye mask.

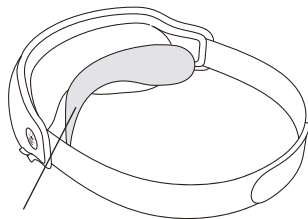


Gel Ice Pack

Caution: Automatic shutdown will happen after 18 minutes in each mode. If you need to stop the massage, press and hold the power button for 1 second to turn off the mask. There will be a long vibration to indicate.

HOW TO CHARGE

- 1 The power box battery indicator level is divided into four compartments. The lights will turn on when the power box is magnetically attached to the eye mask, different light state represents to different power range. When the battery power is too low, the light will keep flashing.

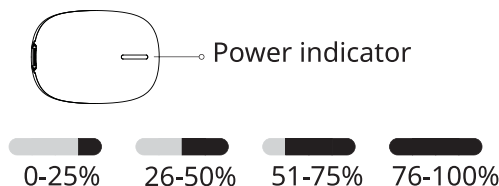

- 2 Before charging, please remove the power box from eye mask and charge it using the USB charging cable. Charging time is approximately 45 min. While charging, the blue indicator light will be flashing, and the light will stop flashing after fully charged. After charging finished, please unplug the charging cable.

TROUBLESHOOTING

Failure	Possible Cause	Solution
Power indicator light is not on	Poor contact between battery and mask	Connect the battery well
	Low battery	Charge the battery
Not working	Low battery	Charge the battery
	Automatically turns off after 18 min of use	Re-start it
Cannot be charged	Poor contact between USB charging wire and power box	Connect the charging wire well
	Your adapter not working well	Change another adapter

WARNINGS

To avoid any damage to person, please follow user instructions details outlined below.

Consult your doctor before use if you are under medical treatment or if you have any of the following conditions:

1. Pregnant women.
2. People who have a body temperature of 38°C or higher.
3. People who have lost sense of their skin temperature.
4. People who have suffered from low temperature burns.
5. People with abnormal body sensation.

People with the following conditions should not use this product.

People who are prohibited from exercise by medical advice.

2. Such as: thrombosis, severe aneurysms, acute venous aneurysms, all kinds of dermatitis, all kinds of skin infections. It is prohibited as it may lead to worsening of these conditions.
3. People who cannot express their will or cannot operate the product by themselves.
4. Children under age of 14 should not use this product.
5. If you feel dizziness, vomiting, throbbing or other physical reactions, please stop using it immediately.

CAUTIONS:

1. Do not use near the fireplace, or electric heaters and other warming appliances.
2. Do not use in a humid environment.
3. Do not use the product in the bathtub, shower, washbasin, swimming pool.
4. Do not place heavy objects on the product after use to prevent damage to the product.
5. Do no disassemble the product.
6. Avoid strong sunlight and corrosive gas and liquid erosion
7. Keep this product out of reach of children.
8. The eye mask is NOT washable. Please use a clean damp cloth to wipe the surface of the mask.
9. Do not allow water to enter the product interior.
10. If it is not used for an extended period, clean it up and put it in a storage bag.



REGISTER E-WARRANTY in 4 SIMPLE STEPS

Please register your warranty online



01

Scan QR Code or
log in to
<https://www.mynt.com.my/e-warranty.html>



03

Register Product
Warranty



02

Attach Product
Invoice



04

Receive Confirmation
Email