




User Manual

MYNT SIMPLY80S MASSAGE GUN

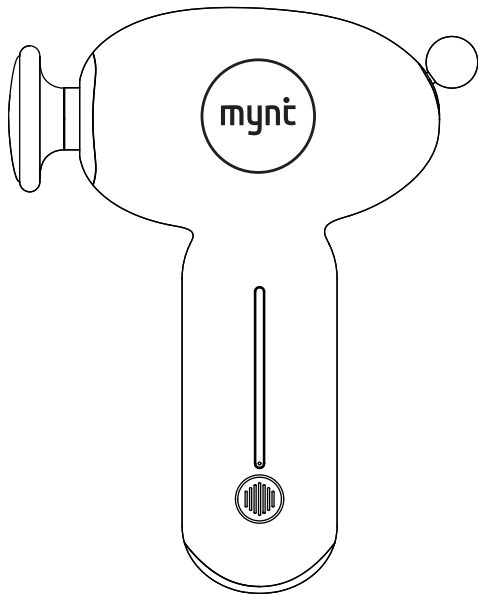

MN-GM01

Please read and keep this manual carefully before using the product.

Thank you for choosing our Massage Gun.

The product has a built-in powerful brushless motor and high-frequency stroke movement, which can perform deep reciprocating beat massage on muscles and other body parts, helping to reduce muscle stiffness and pain, relieve fatigue after exercise, and assist various fitness training. The product is designed to be lightweight, silent, rechargeable, and equipped with a variety of massage heads.

Massage Gun will be your reliable sports partner.

PRODUCT SPECIFICATION


Product Name	Simply80s Massage Gun
Product model	MN-GM01
Power Adapter	5V=2A (Adapter not provided)
Rated Power	20W
Battery Capacity	7.2V/2400mAh/17.28Wh
Charging Time	About 3.5 hours
Auto-Timer	10 Minutes
Penetration Depth	7mm

Our company will continue to improve product performance and appearance, if there is a change, we will not notify separately.

PACKING LIST



Massage Gun

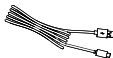

Charging Cable



User Manual



Spherical Head




Cylindrical Head



U-shaped Head

PRODUCT OVERVIEW

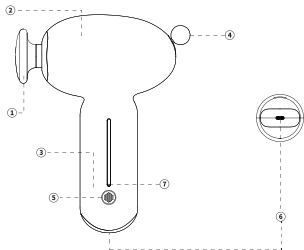

① Massage Head

② Main Body

③ Handle (build-in-battery)

④ Dial Block Ball

⑤ Power Button

⑥ Charging Port

⑦ Indicator Light/
Temperature Indicator Light



Flat Head

For warm compress use
Applicable to limbs and shoulders

Caution: Flat Head is not detachable



Spherical Head

Suitable for impacting deep tissues
such as meridians



Cylindrical Head

Suitable for large muscle groups such
as chest, back, buttocks




U-shaped Head

Suitable for cervical spine and spine

OPERATION GUIDE

STEP 1

Please put massage head on

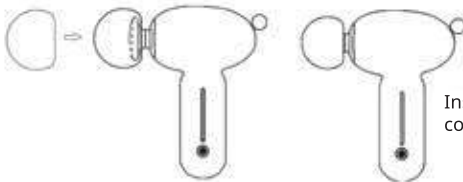

Installation
complete


STEP 2

Press the round power button for 1.5 seconds to turn the device on or off, the white light will illuminate when turned on.


Continue to click the power button to shift gear, which is divided into five gears, including 0 gear for the pause, 1-4 gear gradually increasing in frequency.

Use the Dial block ball to turn on thermal function. The thermal function can only be used after turning on the massager. There are 2 levels of heat, yellow light for 40°C, orange light for 45°C. Use the Dial block ball to adjust the level.

TIPS: Heat therapy best used without massage head.

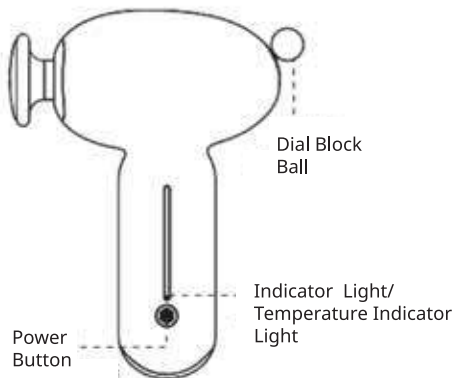

THE LIGHTS TO REMIND

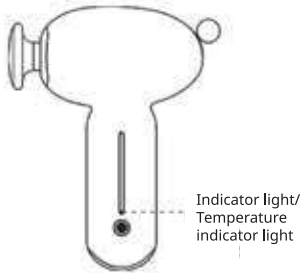

Red Light: Time to charge

White Light: Breathing display (In charging)

White Light always on: Fully charge

Yellow Light: 40°C temperature

Orange Light: 45°C temperature

BATTERY REMINDER

If the light turns red, charge it first;
When charging, the white light shows breathing;
After fully charged, the white light will be always on.

When charging, it cannot be turned on (not broken)

STATEMENT

In order to reduce the risk of malfunction, electric shock, injury, fire, death and equipment damage, Please carefully read all the instructions and warnings in the product manual before use.



Do not touch this part when installing the massage head to start the machine

PROHIBITED

1. It is prohibited for individuals under 18 years of age to use this product. This product is not a toy. Please ensure the safety of children and prevent them from playing with it.
2. If the product accidentally falls into water, it must be repaired before use. Do not use it if not repaired.
3. It is prohibited to disassemble the screws or attempt to modify the product in any way, as it may result in personal or property damage.
4. Using the product on the head, neck, or near the heart is prohibited, as it may pose a life-threatening risk.
5. Using the product directly on the spine is prohibited. Do not use it near the spine without consulting a doctor, as it may cause nerve damage.

6. If you have an injury, it is forbidden to use this product on or near the wound.
7. It is prohibited to use other chargers other than those with a 5V output to charge the product, as it may result in product damage or a fire hazard.
8. It is prohibited for fingers, hair or any other parts of the body to come near the output handheld shaft of the impact shaft, as it may lead to physical injuries.

IMPORTANT MATTERS

If you intend to use the massage gun in any of the following situations, please refrain from doing so without prior consent from a doctor: Pregnancy, diabetic complications (such as neuropathy or retinal damage), wearing a pacemaker, recent surgery, epilepsy or migraines, intervertebral disc herniation, spinal advancement, vertebral detachment or spinal joint stiffness, recent joint replacement surgery or UD, metal implants or any other health related conditions.

Individual with physical weakness should use any vibration equipment under the supervision of nursing staff. These contraindications do not mean that you cannot use vibration equipment or training equipment, but we highly recommend consulting a doctor before doing so.

The current ongoing research is centered around the impact of vibration training on certain medical diseases. This is likely to reduce the categories of contraindications mentioned above. Practical experience has shown that in many cases, incorporating vibration training into a treatment plan is also recommended. However, this should be done under the guidance of a doctor, expert or physical therapist.

WARNINGS AND CAUTIONS

1. This product is not recommended for use by individuals (including children) who are physically impaired, sensory or mentally impaired, or lack of relevant experience and knowledge, unless supervised by someone responsible for their safety.
2. This product is not a medical device and cannot be used for self-treatment. If it is intended to be used to restore muscle damage, please consult your doctor.
3. Do not use it on any hard or bony parts. Please adjust your strength to avoid bruising in the massage area. It is suitable to use the product on the soft tissues of the body without causing any pain or discomfort.
4. Do not use it on wet skin surface, you can gently press and move it over clothes or fabrics.
5. Do not use within one hour before or after eating to avoid adverse effects.
6. Do not use this product when you are drunk or unwell. If you feel unwell during use, please stop using it.
7. Do not let the control panel come into contact with liquids such as water, and do not operate the control switch with wet hands.
8. Do not place any objects in the vents of the product, and do not insert or drop objects into any openings of the product. Blocking the vents may cause the internal components to overheat and leading to malfunction or damage.

9. Do not drop or use the product improperly.
10. When charging, do not pull the power cord or use it to handle the product to avoid wire breakage or poor contact. When unplugging the power plug, hold the plug and pull it out.
11. Using unapproved or incompatible power supplies, chargers or batteries may cause fire, explosion or other dangers.
12. The warm compress function is prohibited for delicate body organs such as the eyes.

COMMON ERRORS AND SOLUTIONS

1. When the machine fails or cannot be started, please check the following items. If the product is still in a faulty state, please send the product to the designated Technical & Servicing Department or authorized technicians for inspection and maintenance.
2. Do not try to dismantle, repair or modify this product by yourself, otherwise it may cause personal injury and property damage.

Failure Phenomenon	Causes of Failures and Treatment methods
Panel light does not light up	<ol style="list-style-type: none">1. Check whether the battery is too low, and long press the power button to see if it can be turned on2. If the battery has been low for an extended period, please charge it for more than 8 hours to check if the battery can be activated. If not, replace the battery with a new one.

<p>The use process stops by itself</p>	<ol style="list-style-type: none"> 1. If the timing is over, please restart the product to continue using it. 2. If the battery is low, please charge it promptly and use it after it's fully charged.
<p>The product runs louder</p>	<ol style="list-style-type: none"> 1. During the running-in period of the new machine, this is normal phenomenon. 2. If there are loose parts inside the product, please send it for repair.
<p>When the product is working, there is a peculiar smell or smoke inside</p>	<ol style="list-style-type: none"> 1. Please stop using immediately, shut down and send it for repair.
<p>The battery life becomes shorter and the battery cannot be charged</p>	<ol style="list-style-type: none"> 1. If the battery is aging, please send it for repair and have the battery replaced. 2. Please send it for repair if the internal parts are faulty.



REGISTER E-WARRANTY in 4 SIMPLE STEPS

Please register your warranty online



01

Scan QR Code or
log in to
<https://www.mynt.com.my/e-warranty.html>



03

Register Product
Warranty



02

Attach Product
Invoice



04

Receive Confirmation
Email